Moisture analyser KERN DBS



Modern moisture analyser with 10 memories and graphic display

| STANDARD | OPTION |
|---|--|
| CAL EXT MEMORY RS 232 USB | GLP PROTOCOL UNIT UNIT UNIT UNIT UNIT UNIT UNIT UNIT |
| KERN | DBS 60-3 |
| Readout [d] | 0,001 g / 0,01 % |
| Weighing range [Max] | 60 g |
| Reproducibility, weight of sample 2 g | 0,15 % |
| Reproducibility, weight of sample 10 g | 0,02 % |
| Display after drying (Display can be switched over at any time) | |
| Moisture [%] = Weight loss (WL) from SW (startweight) | 0 - 100 % |
| Dry mass [%] = Residual weight (RW) from SW | 100 - 0 % |
| ATRO [%] [(SG-RG) : RG] · 100% | 100 - 999 % |
| ATRO [%] (SG : RG) · 100% | 0 - 999 % |
| Weight loss [g] (WL) | Absolute value in [g] |
| Temperature range | 50 - 200 °C in steps up to 1 °C |
| Drying modes | J Standard drying J Drying in levels J Gentle drying J Pre-heat level, can be switched on |
| Switch off criteria | When the set time has expired (1 min - 4 h) When the weight loss per unit of time falls below the target value (30 sec) |
| Recall of measurement | Interval adjustable from 1 sec - 10 min (only in connection with KERN YKB-01N printer or PC) |
| Overall dimensions WxDxH | 202x336x157 mm |
| Net weight | approx. 4,2 kg |
| | |
| Option DKD Calibr. Certificate | 963-127 |





Features

- Drying program
- Previous drying time
- 3 Current temperature
- 4 Unit for displaying the results
- 5 Current moisture content in %
- **6** Drying mode/Status display drying
- Halogen quartz glass heater 400 W
- Observation window over the sample, useful during initial adjustment
- Internal memory for automatic sequence of 10 drying programs and 100 drying processes carried out
- The last value measured remains on the display until it is replaced by a new measurement
- **Sample description** for up to 99 samples, 2 digits, freely programmable, and is printed in the measuring protocol
- Date and time display as standard
- 10 sample plates included
- **Applicationh andbook:** On the internet, you will find a practical application handbook containing many examples, field reports, adjustments and tips for each KERN moisture analyser

Accessories

- Sample plates aluminium, Ø 92 mm. Unit of 80 pieces, KERN MLB-A01
- Roundfi berglass filter e.g. for samples that splash or become encrusted. Unit of 80 pieces, KERN RH-A02
- Temperature calibration set consists of measuring sensor and display device, KERN DBS-A01
- USB interface kit for bi-directional data exchange between moisture analyser and computer. Scope of supply: USB cable, driver, Software Balance Connection, KERN DBS-A02
- Suitable printers see page 130

KERN Pictograms



balance's accuracy with internal adjusting weight (motordriven). Data interface RS-232: To connect the

Internal adjusting (CAL): Quick setting of the



balance to a printer, PC or network Network interface: For connecting the scale to



an Ethernet network. With KERN products you can also use a universal RS-232/LAN converter.



Net-total weighing: weight of tare cup and weight of components memorized in two separate stores.



Weighing with tolerance range: Upper and lower limiting can be programmed individually, e.g. dosing/sorting and portioning. Stainless steel: the balance is protected against



corrosion.

Rechargeable battery pack: rechargeable set.




ACCU

Strain gauges: Electrical resistor on an elastic



deforming body.

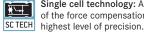

Single cell technology: Advanced version of the force compensation principle with the



Package shipment: The time required to manufacture the product internally is shown in days in the pictogram.



To ensure the high level of precision of your balance, KERN offers the appropriate test weight package for your balance. This consists of the test weight, box and DKD calibration certificate, as proof of its accuracy. The best way to ensure proper balance calibration.

In the extensive KERN test weight range, you will find test weights in the international OIML error limit classes: E1, E2, F1, F2, M1, M2, and M3 with weights from 1 mg to 2000 kg.

The KERN DKD calibration laboratory for electronic balances and weights has been accredited by DKD since 1994 and today is one of the most modern and best-equipped DKD calibration laboratories for balances, test weights and force-measurement in Europe. (DKD = German Calibration Service ~ UKAS)


Your KERN specialist dealer:

Adjusting program (CAL): For quick setting of the balance's accuracy. External adjusting CAL EXT weight required.

Control outputs (optocoupler, _____ digital I/O) to connect relays, signal lamps, SWITCH valves, etc.



GLP/ISO record keeping: of weighing data with date, time and identification-no. Only with printers from KERN.

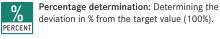

Vibration-free weighing: (Animal weighing Mprogram) Vibrations are filtered out so that MOVE a stable weight is obtained.

Suspended weighing: load support with hook Ē on the underside of the balance. UNDER

Mains adapter: 230V/50Hz in standard version for Germany. On request GB, AUS or 230 V USA version.

Tuning fork principle: A resonating body is $((\mathbf{U}))$ electromagnetically excited, causing it to T-FORK oscillate.




Verification possible: The time required for verification is specified in the pictogram.

Pallet shipment:

The time required to manufacture the product DAYS internally is shown in days in the pictogram.



Memory: Balance contains memories, e.g. for item data, weighing data, tare weights etc. PLU

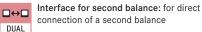



UNIT

Piece counting: Reference quantities selectable. Display can be switched from piece to weight

Weighing units: Can be switched to e.g. nonmetric units at the touch of a key. See balance model. Please refer to KERN's website for more details


Spray and dust protection IPxx: 666 The type of protection is shown by the IP pictogram. For details see the glossary.

| | Battery operation: Ready for battery operation. The battery type |
|------|---|
| BATT | is specified for each device. |

Power supply: integrated in balance. 230V/50Hz in Germany. 230 V More standards e. g. GB, AUS, USA on request.

Electromagnetic force compensation: Coil in N a permanent magnet. For the most accurate FORCE weighings.

DKD calibration possible: The time required DKD for DKD calibration is shown in days in the +3 DAYS pictogram.

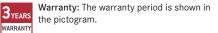

Thanks to the high level of automation, we can carry out DKD calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- DKD calibration of balances with a maximum load of up to 6000 kg
- DKD calibration of weights in the range of 1 mg 500 kg
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DKD calibration certificates in the following languages D, GB, F, I, E, NL, PL

Do you have questions about your scale, the corresponsing test weight or the calibration service ? Your KERN specialist dealer will be pleased to assist you.

KERN - Professional measuring. Measuring technology and testing services from a single source









